

Robust and durable

UDI Code

Latex free



Classic Infusor

Pressure infusion cuff, reusable

/ Classic Infusor

500 ml and 1000 ml Version

With double hook, for fixation of different fluid bag configurations.

- Rapid administration of an infusion solution for compensation of fluid/volume deficiency
- Pressure monitoring and irrigation of the tubing system during invasive, arterial blood pressure monitoring
- Transfusion

3000 / 5000 ml Version

With zipper, allows insertion of different fluid bag configurations.

- Administration of large fluid volumes for irrigation of the operation field

The polyamide material of the Classic Infusor enables easy cleaning & disinfection. Transparent netting for a clear view on the fluid bag and the fluid level. With UDI code for product traceability.

Features




Hand Inflator with Three-Way Stopcock

- Ergonomic hand inflator for inflation
- Three-way stopcock for quick deflation during change of the fluid bag



Manometer

- Pressure range from 0 – 300 mmHg in increments of 10 mmHg
- Robust and durable for accurate pressure monitoring

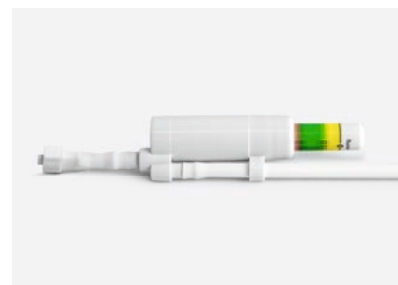

Colour Coded Pressure Indicator

- 360° pressure indicator facilitates pressure monitoring
- Pressure relief valve to prevent over inflation



Order information

Classic Infusor / Reusable

Version	500 ml	1000 ml	3000 / 5000 ml	Box
With hand inflator, three-way stopcock and colour coded pressure indicator	REF 58-02-050	REF 58-02-100	-	1
With hand inflator, three-way stopcock and manometer	REF 58-08-050	REF 58-08-100	REF 58-08-300	1
With Luer Lock connection	-	-	REF 58-01-300	1

Additional information



Application video

This product is manufactured without the use of natural rubber latex, unless otherwise specified.
This product does not contain phthalates which require labelling according to CLP Regulation (EC) 1272/2008.